

# 2200 Series BenchMate

## Vibration-Free Platforms

### Features/Benefits:

- Unmatched price/performance
- Low profile Passive-Air or Active-Air designs
- VibraDamped construction
- Low natural frequency
- Excellent horizontal isolation
- Flatness:  $\pm 0.004$ " over 24" x 24" on -22 and -32 top surfaces
- SPILLPRUF spill management
- Class 100 Cleanroom Compatible – Class 10 available

BenchMate Vibration-Free Platforms enhance the performance of precision tabletop equipment.


The 2210 Series is a Passive-Air design requiring an optional hand pump if an air source is not available. It requires periodic manual leveling.

The 2212 Series is an Active-Air design that requires an air source for use (not included). It is a self-leveling unit and recommended for regular use.

### Ultimate BenchMate

The 2214 Series is very similar to the 2212, however it has the added benefit of multistage horizontal isolation. It is recommended in applications where horizontal vibrations are a problem.

MODEL 2212-01-22

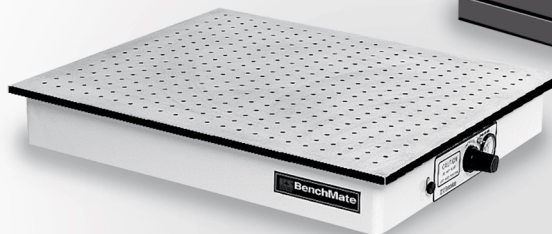

MODEL 2214-01-22



MODEL 2210-04-13




MODEL 2210-01-11




### Performance Data

MODELS 2210 and 2212

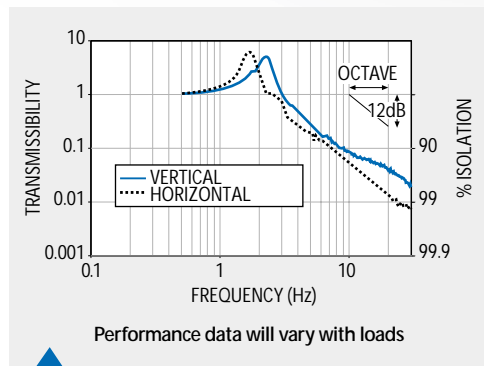

### Performance Data

MODEL 2214

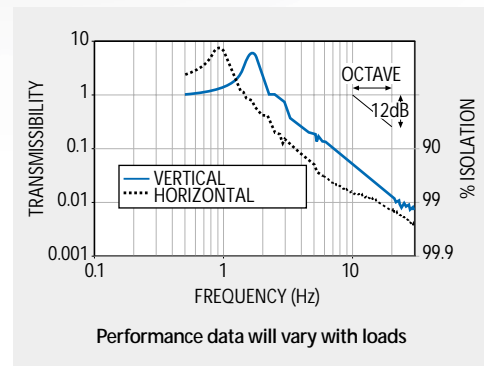

Model 2210 & 2212  
Now Available with Ultraflat  
Breadboard Top  
(-22 & -32 Tops Only)